



**CONFIDENCE. SECURITY. HEM-O-LOK.®**

WECK® Hem-o-lok polymer locking ligation system.

**AT A GLANCE**

- Provides secure ligation, clip locks on to vessels
- Tactile feedback
- Permanent, non-absorbable, inert, non-toxic, radiolucent polymer
- Cool ligation system
- Open or laparoscopic applicators

**CONFIDENCE YOU CAN FEEL**

Confidence. That's what the WECK Hem-o-lok polymer locking ligation system provides. Confidence created by a ligation system designed for clip security—secure in the applicator, locked on to the patient vessel.

**SECURITY BY DESIGN**

Hem-o-lok's unique design and flexible hinge keep the clip firmly seated in applicator jaws. The clip's distinct tactile feedback confirms jaw seating and secure vessel placement.

**COOL LIGATION**


Since Hem-o-lok is a cool ligation system, there's no chance of thermal spread to vital structures, a common concern of many surgeons when they use energy-based ligation solutions.

**SIZE RANGES TO MEET YOUR NEEDS**

Hem-o-lok is available in a range of clip sizes (Medium, Medium/Large, Large, and XL) to ligate up to 10 mm vessels through a 5 mm port or up to 16 mm vessels through a 10 mm port. Our broad range of applicator lengths (including 45 cm endo applicators for extended reach) allows surgeons to deliver clips with control and confidence.




## VESSEL/TISSUE BUNDLE SIZE RANGE

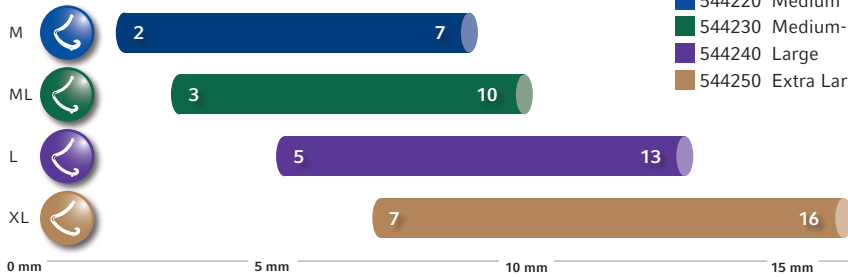

## HEM-O-LOK POLYMER LIGATING CLIPS

- 544220 Medium
- 544230 Medium-Large
- 544240 Large
- 544250 Extra Large

## HEM-O-LOK BEHIND THE CLIP

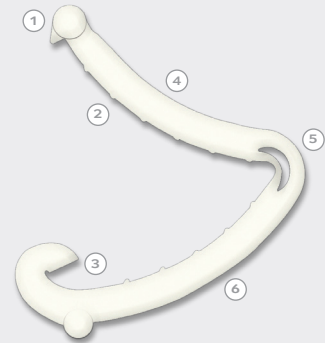

- ① **Bosses** are designed to retain clip in applicator jaws
- ② **Integrated** teeth interface with the vessel and are designed to prevent slippage
- ③ **Locking** mechanism provides tactile feedback and secure closure
- ④ **Bow-shape** allows removal with appropriate instrument
- ⑤ **Hinge** allows flexibility in clip placement prior to clip locking
- ⑥ **Non-Absorbable** polymer is inert, non-conductive, radiolucent, and does not interfere with CT, MRI, or X-ray diagnostics

CLIPS	OPEN APPLIERS/REMOVERS	ENDOSCOPIC APPLIERS/ REMOVERS	AUTOMATIC APPLIERS
 <p><b>M</b></p> <p><b>544220</b> Hem-o-lok® Medium Polymer Clips 6 clips/cartridge 14 cartridges per box, sterile</p>	<p>544113 8" Curved</p> <p>544114 11" Right-Angle (70°)</p> <p>544115 11" Curved</p> <p>544123 11" Remover (for M, ML, L)</p>		
 <p><b>ML</b></p> <p><b>544230</b> Hem-o-lok® Medium-Large Polymer Clips 6 clips/cartridge 14 cartridges per box, sterile</p>	<p>544170 8" Curved</p> <p>544171 11" Curved</p> <p>544172 11" Right-Angle (70°)</p> <p>544123 11" Remover (for M, ML, L)</p>	<p><b>Endo5® ML, 5 mm diameter</b></p> <p>544965 32.5 cm applicator</p> <p>544965L 45 cm applicator</p> <p>544965T 34 cm applicator, take-apart</p> <p>544121 32 cm remover (for ML and L)</p> <p>544121L 45 cm remover (for ML and L)</p> <p>544121T 32 cm remover, take-apart (for ML and L)</p>	<p><b>543965 Auto Endo5®</b></p> <p>15 clips per applicator, 3 applicators per box, sterile, disposable for endo surgery</p>
 <p><b>L</b></p> <p><b>544240</b> Hem-o-lok® Large Polymer Clips 6 clips/cartridge 14 cartridges per box, sterile</p>	<p>544179 11" Right-Angle (70°)</p> <p>544180 8" Curved</p> <p>544181 11" Curved</p> <p>544123 11" Remover (for M, ML, L)</p> <p>544124 11" Remover (for L, XL)</p>	<p><b>Endo10™ L</b></p> <p>544995 32 cm applicator, 10 mm</p> <p>544995L 45 cm applicator, 10 mm</p> <p>544995T 34 cm applicator, 10 mm, take-apart</p> <p>544121 32 cm remover, 5 mm (for ML and L)</p> <p>544121L 45 cm remover, 5 mm (for ML and L)</p> <p>544121T 32 cm remover, 5 mm, take-apart (for ML and L)</p> <p>544130 32 cm remover, 10 mm (for L and XL)</p> <p>544130L 45 cm remover, 10 mm (for L and XL)</p> <p>544130T 33 cm remover, 10 mm take-apart (for L and XL)</p>	
 <p><b>XL</b></p> <p><b>544250</b> Hem-o-lok® Extra-Large Polymer Clips 6 clips/cartridge 14 cartridges per box, sterile</p>	<p>544191 11" Curved</p> <p>544192 11" Right-Angle (70°)</p> <p>544124 11" Remover (for L, XL)</p>	<p><b>Endo10™ XL, 10 mm diameter</b></p> <p>544990 32 cm applicator</p> <p>544990L 45 cm applicator</p> <p>544990T 34 cm applicator, take-apart</p> <p>544130 32 cm remover (for L and XL)</p> <p>544130L 45 cm remover (for L and XL)</p> <p>544130T 33 cm remover, take-apart (for L and XL)</p>	
<p><b>STERILITY STORAGE TRAY</b></p>		<p>544500 capacity: 3 endoscopic 32 cm applicators</p> <p>P37319 capacity: 2-3 endoscopic 45 cm applicators</p>	

TO PLACE AN ORDER CONTACT CUSTOMER SERVICE AT 1.866.246.6990

Teleflex, Endo5, Endo10, Hem-o-lok and Weck are registered trademarks of Teleflex Incorporated or its affiliates. ©2012 Teleflex Incorporated. All rights reserved. 2012-1203