

OLYMPUS

Your Vision, Our Future

VISERA ELITE II

HD 3CMOS AUTOCLAVABLE CAMERA HEAD

CH-S200-XZ-EA

Compact Size and Light Weight



CH-S200-XZ-EA

3CMOS Sensor

- Broad covering of color gamut and enriched color reproducibility in full HD

Special Light Observation

- Select NBI and two modes of IR observation using only one remote button

Small, Compact Designs and Light Weight

- Offers an easy fit in the hand and maneuverability

2× Optical “Power” Zoom

- Magnified observation without image quality deterioration

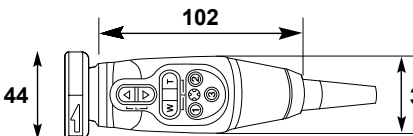
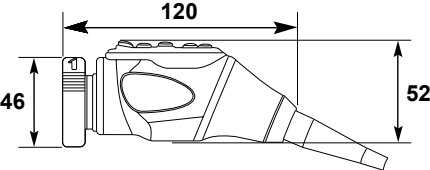
Fiber Mode

- Prevents the moiré effect caused by the combination of camera head and flexible or semirigid scope

Autoclave Compatibility

- Reduced costs compared with other sterilization methods

Specifications

Size	Camera head dimensions (unit: mm)		
	Camera head weight	295 g	
	Cable	ø 6.8 mm × 3.04 m	
Observation	Pickup system	CMOS image sensor (x 3)	
	Magnification ratio	Focal length f = 15.9 to 31.3 mm	
NBI Observation Mode*		Available	
IR Telescopes Observation Mode*		Available	
Electronic Shutter Function*		Available	
Electronic Zoom Function*		Available	
Cleaning/Disinfection/Sterilization	Cleaning/disinfection	Immersible in disinfectant solution	
	Sterilization	Stream sterilization (autoclaving) STERRAD 100S (Short cycle) 100NX (Duo cycle) NX (Standard cycle) V-PRO maX (Non Lumen cycle)	
Classification (Electromedical Equipment)	Type of protection against electric shock	TYPE BF	
	Degree of protection against explosion	The camera head should be kept away from flammable gases	

*For more details, refer to the instruction manual for the video system center.

Specifications, design, and accessories are subject to change without any notice or obligation on the part of the manufacturer.

OLYMPUS

OLYMPUS EUROPA SE & CO. KG
Postbox 10 49 08, 20034 Hamburg, Germany
Wendenstrasse 14-18, 20097 Hamburg, Germany
Phone: +49 40 23773-0, Fax: +49 40 233765
www.olympus-europa.com

E0492496EN · 1.000 · 05/17 · PR · ED