

*Excellent flowrate  
even with inserted instrument*


## **OptiFLOW™**

Continuous highflow, single piece sheath for use with 2.7 mm telescopes • Increased flow dynamics • New distal tip oblique design • Click-Lock • Ideal for the use with ESSURE® device

# OptiFLOW™

## Single Piece, Continuous Flow, Sheath for Hysteroscopy in the Office or Operating room

- Single sheath with Click-Lock for ease of use and added strength
- For use with 30° or 12° R. WOLF 8974 PANOVIEW hysteroscopes
- Increased flow dynamics compared to the traditional dual sheath system
- 45° angle inflow and outflow ports are designed ergonomically for placement of tubing and the most access to the entire working length
- Low profile working channel for a straight approach when placing forceps or the Essure® device
- Oblique tip design is atraumatic for easy access even with a stenotic cervix
- Ports can be removed for cleaning, replacement and sterilization
- One single piece. Ergonomic, easy to use.

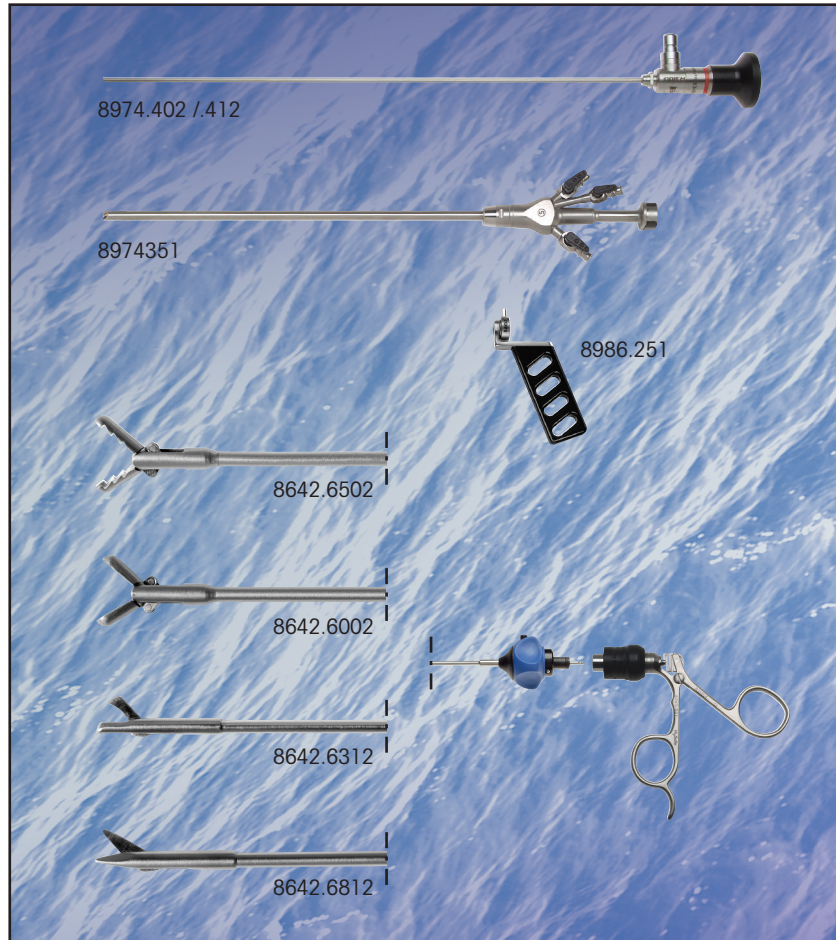

### PANOVIEW Telescope

ø 2.7 mm, working length 300 mm  
 direction of view 30° .....8974.402  
 12° .....8974.412

### OptiFLOW™ Single Sheath

ø 5.5 mm, continuous-flow operating sheath  
 with channel for 5 Fr. instruments  
 for hysteroscopy .....8974351

Ergonomic Handle .....8986.251

### Fibre Light Cable 2.3 m

fibre bundle ø 2.5 mm .....8061.253

### Forceps and scissors:

Micro modular system **HySafe**, rotatable,  
 ø 5 Fr., working length 340 mm

### Grasping forceps

with alligator jaws .....8642.6502

### Biopsy forceps,

spoon .....8642.6002

Biopsy punch .....8642.6312

Micro scissors .....8642.6812

### RIWO-System Tray

400 x 60 x 200 mm (W x H x D)  
 incl. cover and  
 peg mats .....38201.486

