

ENDOCAM Flex HD



At a fraction of the size of traditional endoscopic imaging systems, the ENDOCAM Flex HD significantly reduces your equipment footprint.



Genuine Plug and Play:
Just one plug—connect and go.

Ordering Information

ENDOCAM Flex HD Controller Bundle

ENDOCAM Flex HD Controller (5521101), HDMI/DVI-D cable (103843), power cable (N710171) 55211011

ENDOLIGHT Flex LED Bundle

ENDOLIGHT Flex LED light source (5166001) and power cable (N710171) 51660011

ENDOCAM Flex HD 1080p Camera Head 5521902

RIWO Zoom Lens, with C-Mount thread and snap-on lock, non-autoclavable, f = 13-29 mm, optimized for 1/3" image sensors..... 5261504



spirit of excellence


Be flexible



spirit of excellence

Specifications subject to change without notice.

1430-06-01-0117USA

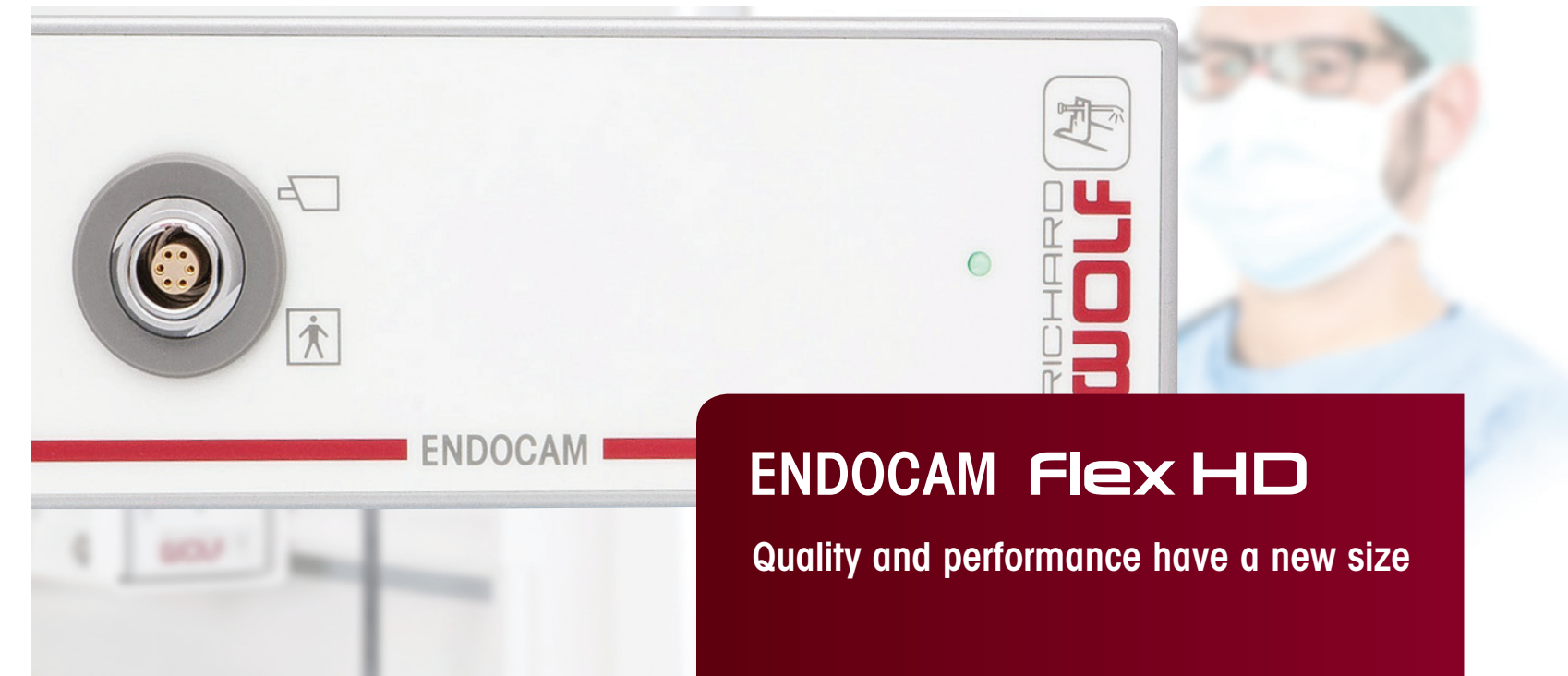

ENDOCAM Flex HD
Quality and performance have a new size

www.richardwolfusa.com

ENDOCAM Flex HD

Be flexible



spirit of excellence



An economical and innovative imaging system, the **ENDOCAM Flex HD** is compact, portable and easy-to-use. This endoscopic camera system offers a high-quality image suited for outpatient and office procedures in areas such as: urology, gynecology, pulmonology and otorhinolaryngology.

With its small footprint, the **ENDOCAM Flex HD** imaging system saves space while delivering brilliant endoscopic images—an optimal solution to meet your needs.

Be flexible

ENDOCAM Flex HD Controller

- Excellent Image Quality—true HD image
- Compact—low-profile controller (150 mm W x 58 mm H x 230 mm D)
- Plug and Play—simple and intuitive setup for easy operation
- Versatile—compatible with conventional rigid and flexible endoscopes



ENDOLIGHT Flex LED Light Source

- Bright—high light output comparable to a 180 W xenon light source
- Environmentally Friendly—power-saving LED technology
- Efficient—up to 30,000 hours of bulb life
- Compact—low-profile light source (150 mm W x 58 mm H x 245 mm D)



ENDOCAM Flex HD Camera Head

- Ergonomic—lightweight (0.28 lbs) 1080p HD camera head
- Snap-on Lock Coupler—13-29 mm parfocal zoom lens
- Two Preset Camera Buttons—white balance (white button) and digital brightness (black button)

The **ENDOCAM Flex HD 1080p system** combines the latest LED technology, intuitive controls and an HD camera head with zoom coupler for optimal image display.